

Opercat®

Self-adhesive dressing pouch

Opercat®

Self-adhesive dressing pouch for the protection of catheters



Product Range (Latex Free)

REFERENCE	POUCH DIMENSIONS	UNIT PER BOX	BOXES PER CARTON	SALES UNIT
0037407	5 x 5 cm	25	13	325
0037417	8 x 7,5 cm	25	13	325
0037427	16 x 5 cm	25	10	250

Fields of Application

DIMENSIONS:	MOST FREQUENT USE:	UNITS:
5 x 5 cm	Catheters for heparinized tracts, normal and Non Return.	Hospitalisation
8 x 7,5 cm	Catheters for yugular and femoral curved, yeyunostomy and gastronomy nutritional drains.	Haemodialysis, Digestive
16 x 5 cm	Catheters like Percath, Hikman, Kerr, yugular and femoral straight. Temporary, lengthy and permanent catheterisations.	Haemodialysis, Surgery, Onco-Hematology, Transplants

- Assures the fixation of catheters
- Quick and easy application and removal process
- Improves the quality of life during catheterisation
- Guarantees the aseptic conditions and sterility of the area



MANUFACTURED BY:
IBERHOSPITEX, S.A.
 Avda. Catalunya, nº 4
 08185 - LLIÇA DE VALL
 (Barcelona) - SPAIN
 Tel.: (34) 93 843 60 34
 Fax: (34) 93 843 61 11
 info@iberhospitex.com
 www.iberhospitex.com




DISTRIBUTED BY:



Opercat® : the best solution

OPERCAT® is a sterile disposable self-adhesive dressing pouch, used for the protection of all types of temporary and permanent catheters.

OPERCAT® is very easy to manipulate by the nursing staff, completely atraumatic for the patient during placement and removal, which improves the asepsis of the area and the patient's quality of life.

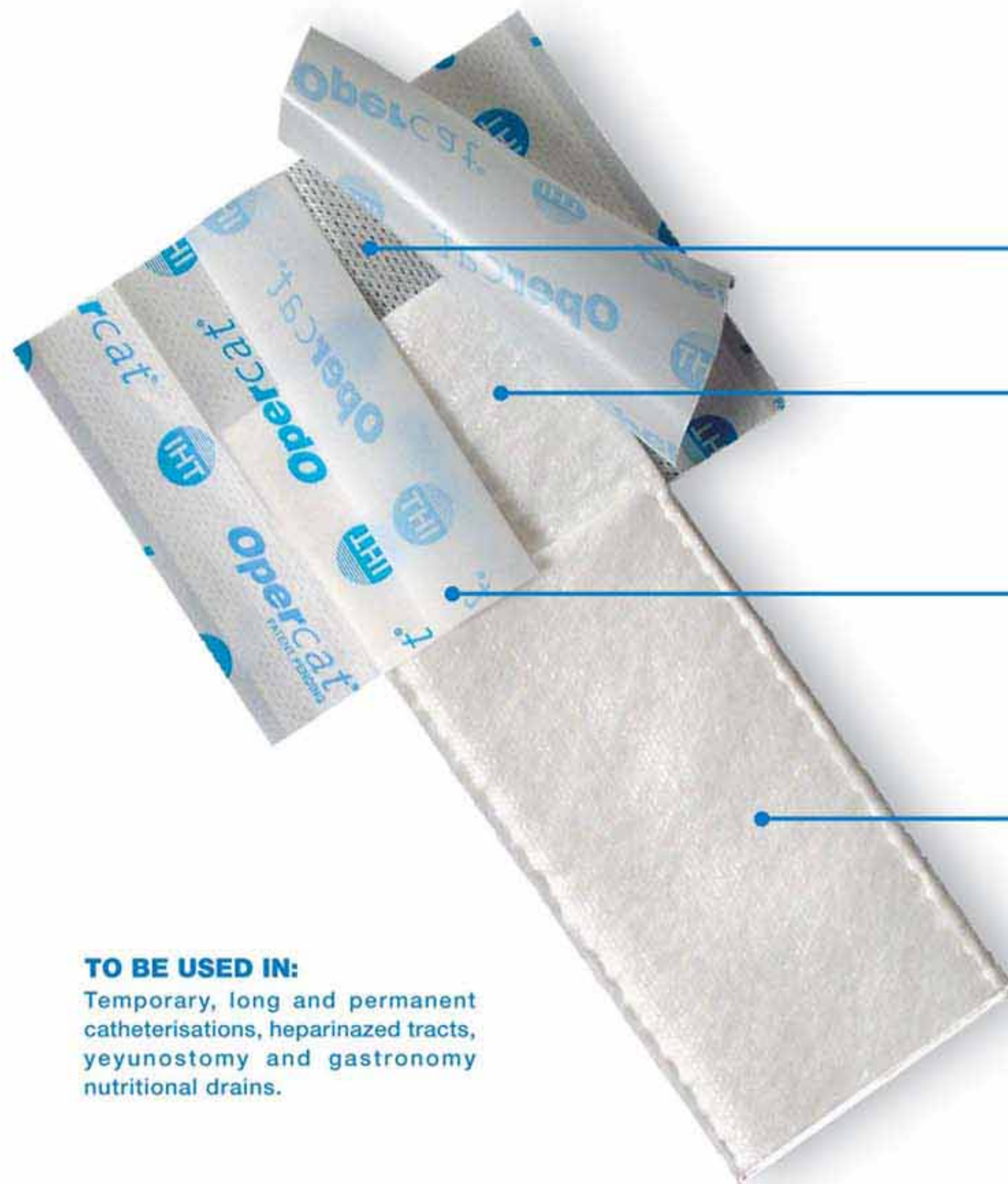

Extensible non woven material with high porosity. Its adhesive power ensures a perfect fixation.

Extension of the same padded material of the pouch (permeable and non adherent) to cover the exit hole of the catheter.

Folded protective sheets to facilitate an easy placement under sterile conditions, protecting the area of the dressing which covers the exit hole of the catheter.

Perspirant pouch with high absorption power, which ensures a perfect padding of the catheter and acts as a barrier for the microorganisms. Completely covered with a permeable and non adherent inner and outer coating.


TO BE USED IN:

Temporary, long and permanent catheterisations, heparinized tracts, yeyunostomy and gastronomy nutritional drains.

APPLICATION

After the placement of the catheter inside the pouch and the removal of the protection sheets, the non-woven adhesive strip has to be fixed on the patient's skin, achieving a double function:

- To fix the pouch to the patient's body
- To keep the catheter secure under complete asepsis

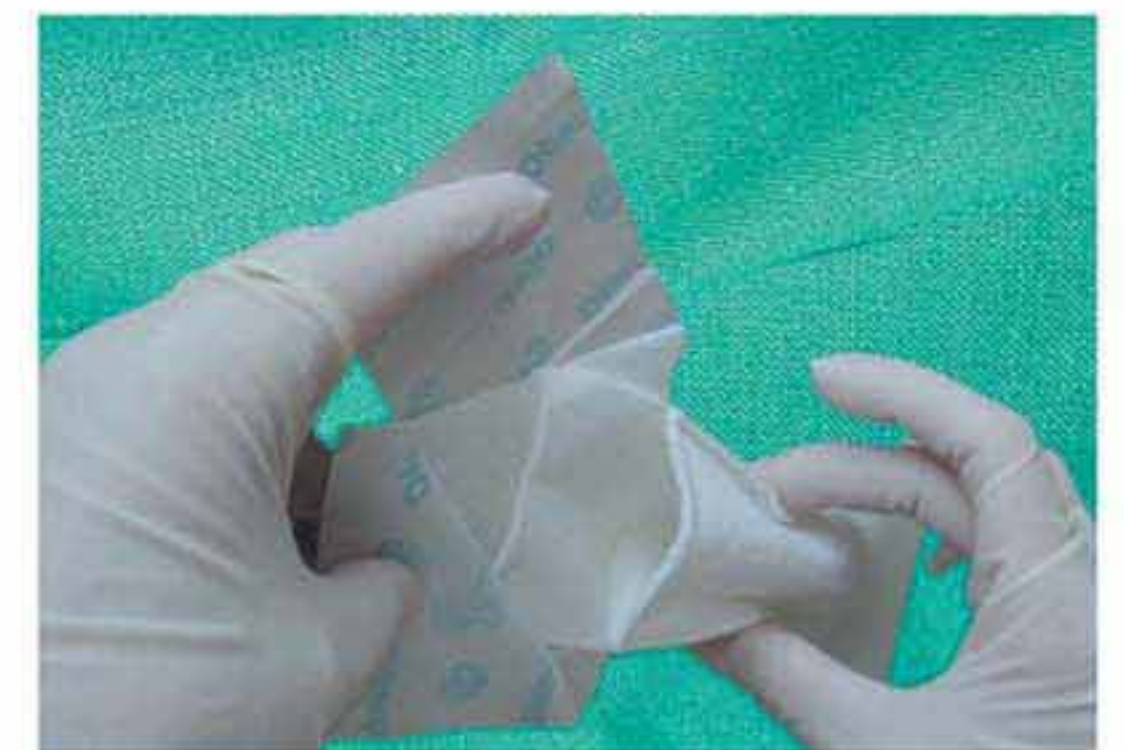

Quick and easy placement, ensuring a secure fixation and protection of the catheter.



Makes the patient feel comfortable during the catheterisation period, and facilitates the cleanliness even at home.



No displacement of the catheters. Easy and comfortable removal of the dressing without scissors or instruments.